


NO LATEX

CODE : **PT PVC 130SL**
FM code : **LIB55**

Extension line 130 mm with T

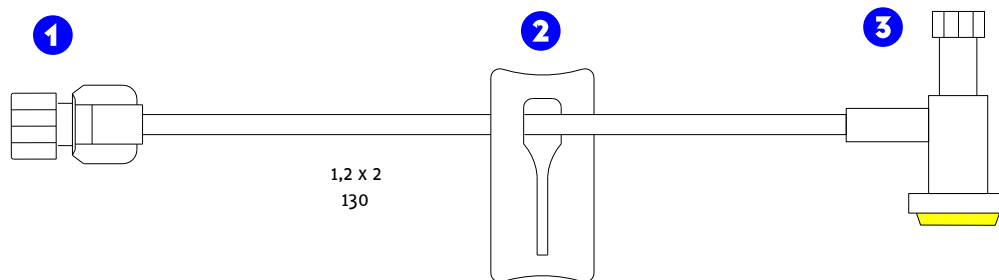

- 1** DEHP-FREE PVC **Female luer lock connector** with vented cap made of PP.
- 2** White PP **sliding clamp**.
- 3** **Luer Male T connector** thin Latex-free injection site made of ABS with cap made of PE.

Tubing is made of soft PVC.
(Tubing Lengths expressed in mm.)

NO LATEX



EffeEmme
FABBRICAZIONI MEDICALI

CODE : **PT PVC 130SL**
FM code : **LIB55**

Extension line 130 mm with T



CE risk class : IIa

Certification Body:

F.M. S.p.A. is TÜV SÜD certified, both for the CE (CE0123) mark, and the UNI EN ISO 9001/ISO 13485 standard.

Environment and Production Process:

The assembly is made in a 100.000 class clean room. Each material is tested for biocompatibility according to ISO 10993 standard.

Packaging:

Blister made of 2 different materials: plastic film and medical - grade paper.

Box size: cm. 60 x 31,5 x 28 h.

Packaging Unit: 500 pcs. / box (10 small boxes of 50 pcs. each).

Easy and safe opening of boxes, through a tab, to avoid the usage of cutters.

Storage and Use:

The circuit sterility is assured in a free from breakage package. To be used immediately after having removed the caps. The expiry date refers to the product in unopened and properly stored package. Recommended storage temperature: from 0 ° to + 40 °C.

Sterilization:

Sterile by EO.

Process according to UNI EN ISO 11135, EN 556. Sterility guaranteed for 5 years. Not to be resterilized.

Product Shelf Life: 5 years from date of manufacture.

Compliance:

Standard applicable to the product: ISO 13485, ISO 10993, ISO 11135, EN 20594, ISO 11607, EN 1041, EN ISO 14971.

Produced by:

F.M. S.p.A.
Via Farini, 65 bis
13043 CIGLIANO (VC)
ITALY

EffeEmme
FABBRICAZIONI MEDICALI