


DEHP FREE
NO LATEX



Code	FM code	Length mm
ML PDF900MF	LIB60	900
ML PDF1200MF	LIB56	1200
ML PDF1500MF	LIB57	1500
ML PDF1800MF	LIB58	1800
ML PDF2100MF	LIB59	2100

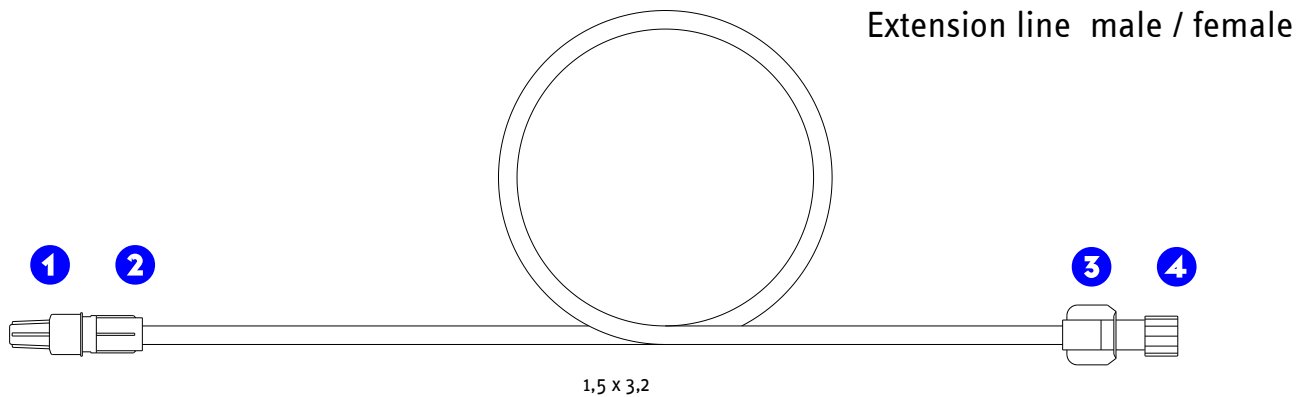

Extension line male / female

- 1** Vented cap made of PE.
- 2** Male luer lock connector made of ABS.
- 3** Female luer lock connector made of DEHP FREE PVC.
- 4** Cap made of PE.

Tubing is made of DEHP-FREE PVC 86 ShA.
(Tubing Lengths expressed in mm.)

CE risk class : IIa

Certification Body:

F.M. S.p.A. is TÜV SÜD certified, both for the CE (CE0123) mark, and the UNI EN ISO 9001/ISO 13485 standard.

Environment and Production Process:

The assembly is made in a 100.000 class clean room. Each material is tested for biocompatibility according to ISO 10993 standard.

Packaging:

Blister made of 2 different materials: plastic film and medical - grade paper

Box size: cm. 58 x 22 x 26,5 h.

Packaging Unit: 100 pcs. / box

Easy and safe opening of boxes, through a tab, to avoid the usage of cutters.

Storage and Use:

The circuit sterility is assured in a free from breakage package. To be used immediately after having removed the caps. The expiry date refers to the product in unopened and properly stored package. Recommended storage temperature: from 0 ° to + 40 °C.

Sterilization:

Sterile by EO.

Process according to UNI EN ISO 11135, EN 556. Sterility guaranteed for 5 years. Not to be resterilized.

Product Shelf Life: 5 years from date of manufacture.

Compliance:

Standard applicable to the product: ISO 13485, ISO 10993, ISO 11135, EN 20594, ISO 11607, EN 1041, EN ISO 14971.

Produced by:

F.M. S.p.A.
Via Farini, 65 bis
13043 CIGLIANO (VC)
ITALY